



Auto Seal has arrived.

EasyLife is so extraordinary, it practically fits itself.

Sometimes it's okay to expect immediate results. That's why we developed EasyLife with revolutionary Auto Seal technology.

We named it EasyLife. You'll call it your favorite set-up mask. With its patented Auto Seal technology, the EasyLife nasal mask nearly eliminates the need for manual mask adjustments. Its innovative, lightweight design features a unique dual-cushion construction – the inner cushion creates an instant, self-adjusting seal, while the outer cushion provides comfortable support. The only manual steps necessary are intuitive headgear adjustments, making fitting as easy as possible. Helping you fit 100% of your patients, 100% of the time is, quite simply, easier with EasyLife.

Hitwin
HEALTHCARE

Respiratory • Sleep Therapy • Cardiology

PHILIPS

RESPIRONICS


Forehead pad adjusts automatically – eliminating a traditional fitting step – and requires minimal headgear force.



Only four parts:
mask frame, headgear,
outer support cushion,
inner seal cushion.

Angled exhalation micro
ports make operation
quieter and redirect air
away from a bed partner.

Outer support cushion lets the mask
rest lightly and comfortably against
the patient's face, eliminating the
need to over-tighten the headgear.

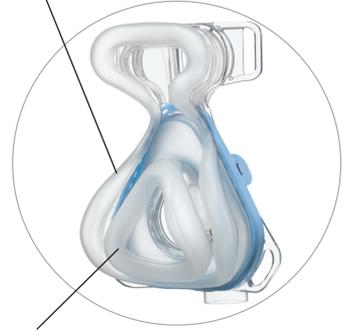

Inner seal cushion creates an instant,
self-adjusting seal.

Ordering information					
EasyLife mask	Small	Medium	Medium wide	Large	U.S. HCPCS
Mask with headgear	1050001	1050002	1050003	1050004	A7034, A7035
Mask only	1050011	1050012	1050013	1050014	A7034
DuoPack (includes headgear and two cushion sets* of the same size)	1050006	1050007	1050008	1050009	A7034, A7035, A7032
FitPack: includes headgear and small and medium cushion sets*		1050078	N/A	N/A	A7034, A7035
FitPack: includes headgear and medium wide and large cushion sets*	N/A	N/A		1050079	A7034, A7035
Replacement parts				Part number	U.S. HCPCS
Small cushion set*				1050091	A7032
Medium cushion set*				1050092	A7032
Medium wide cushion set*				1050093	A7032
Large cushion set*				1050094	A7032
EasyLife headgear				1050095	A7035
EasyLife sizing gauge				1050096	N/A

*A cushion set includes one inner seal cushion and one outer support cushion.
Respironics, Auto Seal, DuoPack, EasyLife, and FitPack are trademarks of Respironics, Inc. and its affiliates, all rights reserved.

Please visit www.philips.com/easylife



LEADING DISTRIBUTOR

Hitwin Healthcare Pvt. Ltd.
7/10, Jaganathan Nagar 2nd Main Road,
Arumbakkam, Chennai, India - 600106.
Ph: 91 44 4684 5555 (20 Lines), 91 44 4233 3313
Mob: 91 93441 66778

CHAT WITH US

